

## lifetime limited warranty

Your Drive branded product is warrantied to be free of defects in materials and workmanship for the lifetime of the product for the original consumer purchaser.

This device was built to exacting standards and carefully inspected prior to shipment. This Lifetime Limited Warranty is an expression of our confidence in the materials and workmanship of our products and our assurance to the consumer of years of dependable service.

This warranty does not cover device failure due to owner misuse or negligence, or normal wear and tear. The warranty does not extend to non-durable components, such as rubber accessories, casters, and grips, which are subject to normal wear and need periodic replacement.

If you have a question about your Drive device or this warranty, please contact an authorized Drive dealer.



© 2014 Medical Depot, Inc. All rights reserved.  
Drive is a trademark of Medical Depot, Inc.  
Port Washington N.Y. 11050 USA Made in China

## deluxe steel drop arm commode

item # 11125KD-1 & 11125KD-4



EU Authorized Representative  
Drive Medical LTD  
Ainley's Industrial Estate  
Elland, West Yorkshire,  
United Kingdom HX5 9JP

[www.drivemedical.com](http://www.drivemedical.com)


A large, vertical version of the Drive logo is positioned on the left side of the page. The word "drive" is written in white, lowercase, sans-serif font against a dark blue background. A small yellow square is positioned above the letter "i". A trademark symbol (TM) is located to the upper right of the word.

[www.drivemedical.com](http://www.drivemedical.com)


## commode assembly

1. Remove contents from box. Depress brass push pins and insert bucket frame section and frame support bar into arm section as shown below. The push pins will snap into place when frame is securely in place.

Note: Make sure the back of the bucket frame (the side with the back rest receptacles) corresponds to the back of the armrest section.

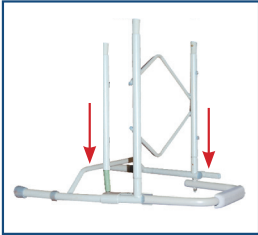

2. Repeat for second arm section.

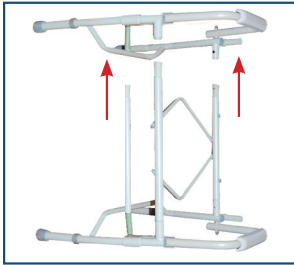

3. Insert backrest into backrest receptacles in commode bucket frame. The brass push pins will snap into place when secure.



4. Align the toilet seat clamps onto the back frame rail. (Both clamps should sit outside the pail brackets). Push down on seat until clamps snap firmly into place. Repeat for the lid. To remove, grasp the back of the seat and lift.



5. Install commode bucket by dropping into bucket frame. Reverse to remove.



## adjusting leg height

1. Before use, adjust height of each leg to your specific needs, making sure all legs are at the same height. When desired height is reached, lock legs into place by aligning buttons into proper holes
2. To be sure legs are locked into place by pushing down on or by pulling legs.

## additional precautions

1. Toilet seat must be in down position before use. Failure to do so may result in breakage of commode and possible injury.
2. Do NOT install commode without first reading and understanding this assembly and instruction sheet. If you can not understand these instructions, contact a healthcare, dealer or technical professional before attempting to install. Otherwise, injury or damage may occur.
3. Inspect all parts for shipping damage before assembly. If there is shipping damage - DO NOT USE. Contact dealer/carrier for further instructions.
4. Make certain that the height adjusting snap buttons fully protrude through the same respective hole of each leg extension. This will ensure that the leg extensions are securely locked in position and an even height is achieved.
5. Ensure ALL screws, nuts and/or bolts are tightened.
6. Check the rubber tips on the leg extensions for rips, wear or if they are missing. Immediately replace any or all if any of these imperfections exist.
7. All four leg extensions with rubber tips MUST touch the floor simultaneously at all times.
8. Users with limited physical strength should be supervised or assisted while using commode.
9. To operate the drop arm, push the lever forward to release and drop the arm. To return the arm rest to the upright position, pull the arm up until the lever latch clicks into place. Caution, do not lean or push down on the arm rest until the latch has been checked to insure the latch is in place and will hold weight.
10. Use the arm supports only for assistance. DO NOT attempt to use the arms to support full body weight.
11. Arms should be unlatched and dropped to assist in transfer onto/off the device